BellCROSS™


CE

VM-EP01

By Visiomed®

A**1**

37,7 x 7,6mm

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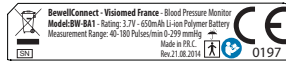COMPLIANT WITH EUROPEAN STANDARDS

CONTENTS

1. NOTES ON SAFETY11

2. FEATURES12

3. USE13

3.1 DESCRIPTION.....13

3.2 PRECAUTIONS BEFORE USE13

3.3 USER GUIDE13

4. HOW TO MAINTAIN AND STORE THE UNIT13

4.1 INSERTING BATTERIES13

4.2 CLEANING THE APPLIANCE14

4.3 MAINTENANCE14

5. HOW TO PREVENT A MALFUNCTION14

6. SPECIFICATIONS.....15

The manufacturer keeps the right to modify without any preliminary opinion technical specifications of the product.

Dear customer,

Thank you for choosing the VM-EP01 BELLCROSS Pedicure. This manual provides important information to help you use the system properly. Before using this product, please read the following contents thoroughly and carefully.

If you have other questions regarding this product, please contact the local customer service or place of purchase.

1. NOTES ON SAFETY

- Use on dry skin only (not wet).
- Do not apply to the same area for more than 3-4 seconds at a time.
- Discontinue use immediately in the event of irritation or discomfort.
- Do not use on chapped, sore or irritated skin.
- Do not immerse the device in water.
- Do not use this device if suffering from diabetes or poor circulation.
- Keep long hair and clothing away from the device during use.
- This appliance must always be placed in a clean and dry place. Do not use this appliance near bathtubs showers, basins or other vessels containing water. Keep it dry.
- Designed for smoothing coarse skin on the feet and heels. Do not use this device anywhere else on the body.

BellCROSS™ EN

USER MANUAL

VM-EP01

By **Visiomed**

- Do not use on broken or bleeding skin.
- Do not place in wet or high temperature.
- Keep out of reach of children.
- Do not disassemble the body.
- Use this appliance only for the use for which it is designed as described in this manual.
- This appliance may be used for personal use at home. Do not use the unit for other purpose.
- Use this appliance in an ambient temperature range of 10 to 40°C.
- Do not expose this appliance to extreme temperature conditions > 60°C ou < -20°C.
- Do not use this appliance at a relative humidity of more than 85%.
- Do not expose this appliance to sunlight.
- Do not expose this appliance to electric shocks.
- Do not use this appliance outdoors.
- Never drop the appliance.
- Follow the maintenance instructions specified in this manual.
- Do not attempt to open the device. In case of problems, contact your seller.
- Discontinue use of the device in case of anomalies or malfunction.
- This appliance is not designed for use by people (including children) with reduced physical, emotional or mental abilities or those lacking experience or knowledge unless used in the presence of someone responsible for their security, under surveillance or following instructions about how to use the appliance.
- Do not use any battery other than the mentioned batteries. Do not attempt to recharge non-rechargeable batteries. Do not throw them into a fire.
- Remove the batteries if the appliance is not used for a long period or regularly.
- Emergency treatment: Power off.

Thank you for choosing the VM-EP01 BELLCROSS electronic Pedicure

IMPORTANT: ENSURE THE POLARITY OF THE BATTERIES IS CORRECT. INCORRECT POLARITY CAN CAUSE DAMAGE AND COMPROMISE THE WARRANTY OF YOUR APPLIANCE.

2. FEATURES

- Unique rollers spin 360° to gently remove hard skin
- Safely switch-lock, avoid unwanted operated
- Exfoliating rollerX2
- Get Smoothest Feet

3. USE

3.1 DESCRIPTION

View diagram A page 2

- 1 Protective Cover (Remove the protective cover before use)
- 2 Exfoliating Roller
- 3 Release Button
- 4 Switch
- 5 Battery cover

3.2 PRECAUTIONS BEFORE USE

- Do not use on broken or bleeding skin.

3.3 USER GUIDE

- Wash and clean the skin and dry fully.
- Make sure the Roller is placed into the unit.
- Remove the protective cover from the unit and it is ready to use.
- Push switch up to turn unit on.
- Gently move to callus remover back and forth or side to side over the callus/corn/hard skin area for 2/3 seconds at a time. Do not press hard and gently glide over the desired area. The unit will stop if it is pressed top hard against the skin. Stop to check if you have achieved the softness you desire. If not, repeat again 2/3 seconds at a time, continue like this until the level of smoothness desired is obtained.
- Turn the unit off.
- Rinse the skin or use a wet towel to get rid of the excess dead skin particles.
- Reattach the protective cover.

4. HOW TO MAINTAIN AND STORE THE UNIT

4.1 INSERTING BATTERIES

1. Alkaline batteries are recommended
2. For long periods of non-use, please remove battery

Operation of battery case and the power switch

1. Insert batteries into callus remover


2. Tighten the battery cover

3. Turn on the callus remover by pushing the switch up

4.2 CLEANING THE APPLIANCE

- Remove the protective cover from the head.
- Clean with brush. Do not turn the appliance on while brushing.
- Do not submerge the device in water.
- Turn the unit off.
- Wipe the body of the unit with a clean dry cloth.
- Replace the protective cover.

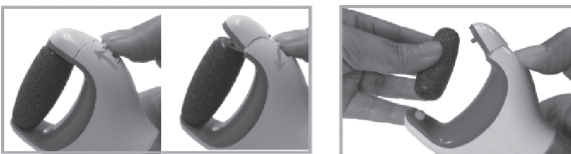4.3 MAINTENANCE

- Do not place the appliance in an area with excessive heat or moisture.
- Keep out of children's reach

5. HOW TO PREVENT A MALFUNCTION

- Avoid dropping and heavy impact.

- Do not put the unit in water.
- Do not attempt to disassemble the unit.
- The appliance might not meet its performance specifications or cause safety hazard if stored or used outside the specified temperature and humidity ranges in specifications.

6. SPECIFICATIONS

Product name	BELLCROSS Pedicure
Model	VM-EP01
Voltage	DC 6V
Power	4 W
Battery	4 x AA (included)
Size of the unit	165 x 72mm
Weight	128 g (not include batteries)
Speed	1700-1900 trm
Material	ABS

Notes : The roller is not under warranty.

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EN: Visiomed* will repair or replace this product free of charge in the case of defective parts or manufacturing defects, in accordance with the conditions mentioned below as follows:

DURATION: 24 MONTHS RETURN TO WORKSHOP

LIMITS AND EXCLUSIONS: This guarantee concerns only the original final purchaser. A purchase invoice, or another proof of purchase, with this guarantee card will be required to obtain an after-sales service, in accordance with this guarantee. This guarantee card will not be extended to another person only the original final purchaser. This guarantee becomes void if the serial numbers on the product are modified, replaced, illegible, absent, or if repair has been carried out by a service not approved, including the user.

This guarantee covers only the defects of the material or parts, occurring during normal use of the product. It does not cover the damage caused during the transport of the apparatus, causes due to repairs being carried out by the distributor, by any modifications undertaken, any connection of equipment not approved by Visiomed*, or causes contrary to those written in the user manual or notice. Moreover, the present guarantee does not cover damage due to falls, bad handling, bad installations, damage by fire, floods, lightning, or any other natural disaster. This guarantee does not cover the packing of the material, the accessories, the defects caused by commercial exposure of the product, show room, sale space, demonstration etc... Normal maintenance, cleaning and the replacement of parts where wear is normal, are not covered by the terms of this guarantee. Visiomed* and its representatives and agents will not in any case be held responsible for any damage and consecutive damages due to the mishandling of this product. This guarantee is the only valid one at Visiomed*, any other guarantee (commercial guarantee) except this one will not be taken into account.

IMPORTANT: During the guarantee period if you are dissatisfied with the repairs of this product, please contact the Visiomed* customer service.

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EN: Disposal of Waste Equipment by Users in Private Households in the European Union.

The symbol on the product or on its packaging indicates that this product must not be disposed of with your other household waste. Instead, it is your responsibility to dispose of your waste equipment by handing it over to a designated collection point for the recycling of waste electrical and electronic equipment. The separate collection and recycling of your waste equipment at the time of disposal will help to conserve natural resources and ensure that it is recycled in a manner that protects human health and the environment. For more information about where you can drop off your waste equipment for recycling, please contact your local city office or your household waste disposal service.

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© 07/2016 BW-EP01
EN: *BewellConnect® is a brand registred by Visiomed Group SA. The constructor keeps the right to modify the photographs and characteristics without previous notice. Registred design. The logos and trademarks are the property of their respective owners.*